


Baylis
MEDICAL

Nykanen RF Wire



Nykanen RF Wire

FOR CREATING A CONTROLLED PUNCTURE

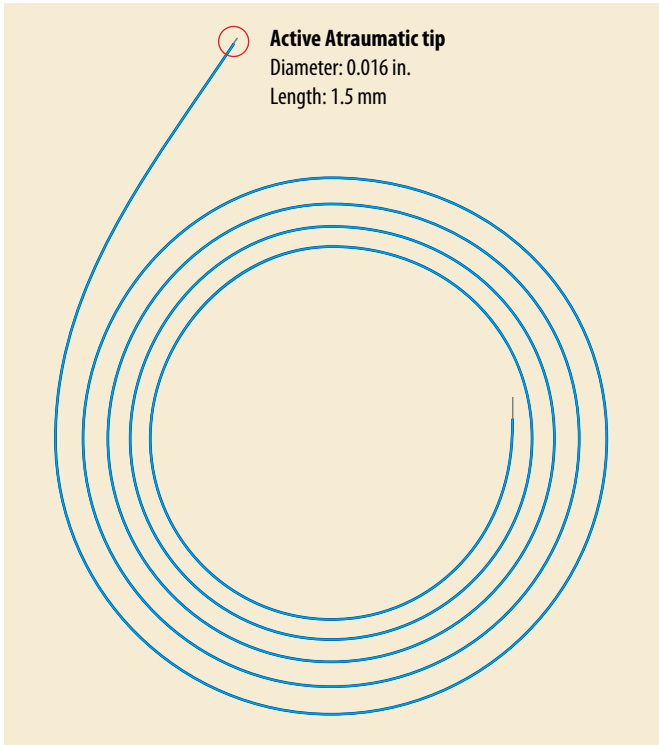


The Nykanen RF Wire is designed to create a controlled puncture by delivering radiofrequency energy. Physicians have used the wire when treating conditions associated with pulmonary atresia, hypoplastic left heart with intact atrial septum, and vascular occlusions.

Excellent Handling

- The wire handles like a coronary guidewire
- The wire has a flexible distal region that can be manually bent and shaped with good memory
- Its rounded, atraumatic tip will not mechanically perforate tissue
- Its PTFE insulation facilitates smooth advancement through guiding catheters

Designed For Quick Exchange

- The proximal end of the wire has no fixed connector, which allows the BMC Coaxial Injectable Catheter to be loaded over top of it
- Its ample length of 265 cm allows for the 145 cm ProTrack™ microcatheter to be loaded onto the Nykanen wire outside of the patient

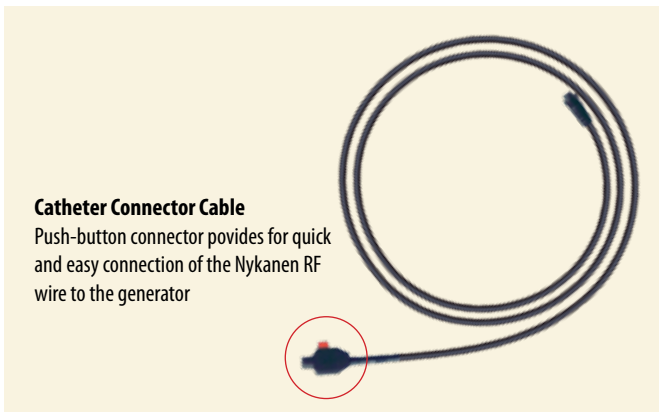

Specifications

Nykanen RF Wire	RFP-265
Body Outer Diameter	0.024 in
Overall Length	265 cm
Catheter Connector Cable	RFP-101
Overall Length	10 ft



© Copyright Baylis Medical Company Inc., 2007. Baylis Medical Company Inc. reserves the right to change specifications or to incorporate design changes without notice and without incurring any obligation relating to equipment previously manufactured or delivered. ProTrack™ and Baylis Medical Company Inc. logo are trademarks or registered trademarks of Baylis Medical Company Inc. in the United States of America and/or other countries.

CAUTION: Federal Law (USA) restricts the use of this device to or by the order of a physician.

Patents Pending and/or issued

CAR1008



Baylis Medical Company Inc.
5959 Trans-Canada Highway
Montreal, QC Canada H4T 1A1

Tel. : (514) 488-9801 / Fax : (514) 488-7209
www.baylismedical.com / info@baylismedical.com